



Bruker SPR Remote Access API



Revision H
Control Software version 3.9.2 and higher
2026-08-18

Overview

The Bruker SPR Control Software's remote API allows client frameworks to connect to the Control Software remotely and run its most important operations like starting batch jobs or maintenance operations.

The implementation is based on the Microsoft Windows Communication Foundation (WCF) technology which allows us to publish and configure several different endpoints to access the API by just changing the Control Software's configuration file. This allows writing clients to use a wide range of different languages, connections, and protocols to access the API. To keep the API as flexible as possible, it does not support callbacks. Updates about the system's status must be obtained using polling.

As the API is specific for the control software, it is highly recommended to have a brief understanding of what you can do with the control software, what the SPR machines of Bruker are capable of, and how they (in principle) work. For that it might be helpful to investigate the manual of the Control Software. The manual is integrated into the software and can be opened by clicking on description in the right-side bar.

Remote-Control Mode

To use the remote API, some parameters in the Control Software's configuration file might need to be adjusted. For SPR 24 and 32 Pro devices the configuration file is *MASS1.ini* and for SPR #64 devices it is *SPR64.ini*.

1. Ensure that `MultipleRunsetMode` has the value `True`.

Example:

```
[Application]
MultipleRunsetMode = True
```

- 2.

- a) For SPR 24 and SPR 32 Pro devices ensure that `DeckFile` has the value

```
..\Racks\DeckMASS1_Orbitor.xml
```

Example:

```
[LiquidHandling]
DeckFile = ..\Racks\DeckMASS1_Orbitor.xml
```

- b) For SPR #64 devices ensure that `DeckFile` has the value

```
..\Racks\DeckSPR64_exchangeable_P2_P4.xml or
```


```
..\Racks\DeckSPR64_exchangeableP4.xml
```

Example:

```
[LiquidHandling]
DeckFile = ..\Racks\DeckSPR64_exchangeableP4.xml
```

Before a client can connect to the API, the Control Software must be switched to the remote-control mode. This mode disables most buttons of the user interface so a user in front of the device cannot interfere with the remote client's commands. The user can only check the status of the device and look at the data. The control software opens a server to which a client can connect.

The current user is logged off from the control software when the server is started since the remote-control mode uses a dedicated user account with special permissions.

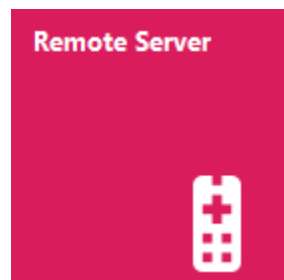

The remote-control mode is started by clicking the “Remote Server” tile on the Control Software’s home screen.

You will get to a new page where you can see a graph with the application’s recorded data and a log table showing messages that are exchanged between server and client. Once the server is started, all navigation buttons of the user interface are blocked, and it is not possible to leave the page while the server is running.

Message boxes that are displayed during the regular operation of the Control Software are intercepted and handled by the server to avoid the software getting blocked while a dialog is shown. Dialogs that allow more than one possible answer can be configured in the file “RemoteServiceDialogResults.xml” which is part of the Control Software’s installation.

Simulation Mode

The purpose of the simulation mode is to allow users to plan their assays without using a real machine. It may also be useful for the developers of lab automation frameworks who integrate the SPR device into their framework. The software is automatically started in simulation mode if no machine is connected.

There is also a feature to simulate errors. It is implemented as batch command (“Throw Error” which you will find in the Method Editor) which can also be used for nonremote testing. The command has no parameters to select. The command will throw an exception when it is executed during a batch job and cause the job to crash with an error message. To use the command, you first must enable it. Go to the folder “<installation directory of the control software>/Racks” and edit the file “CommandPrototypesSPR32.xml”. At the end of this file you will find a “BatchMethodElement” with the “<Name>Throw Error</Name>”. Remove “<!--” before “<BatchMethodElement ...>” and the “-->” after “</BatchMethodElement>” and save the file. Afterwards you must restart the software. To simulate an error, you must create a new method with the Method Editor and add the command “Throw Error” to this newly created method. You can also add additional commands before and after the “Throw Error” command and you can add multiple “Throw error” commands.

Endpoint Configuration

The endpoints to which clients can connect are configured in the Control Software’s *app.config* file. Endpoints can be configured for different transfer protocols (TCP or HTTP). Clients can either run locally on the same PC as the Control Software or access the API via a network connection. The base addresses for the endpoints can also be configured. By default, the base address for HTTP endpoints is 9001, resp. 9002 for TCP connections.

Once the remote-control server of the Control Software has been started you can check that it is running by entering the URI <http://localhost:9001/bruker-spr> into a web browser’s address field on

the computer on which the Control Software is running. You can also do this from a different computer if "localhost" is replaced with the IP address of the Control Software's PC. Of course, port 9001 must be open to the outside to do this.

You can also use the JSON endpoint to transfer commands directly from the browser to the control software. Some examples:

- <http://localhost:9001/bruker-spr/json/GetOperationMode> gets the current operation mode as an integer. The meaning of the value can be picked from the table in the Interface Description section.
- <http://localhost:9001/bruker-spr/json/GetNamesOfMaintenanceProcedures> will return a list of names of all maintenance commands that can be run remotely in JSON format.
- <http://localhost:9001/bruker-spr/json/RunMaintenanceProcedure?procedureName=Reset%20System> starts the maintenance command to reset the device's stepper motors and pumps. This only works if the Control Software is either idle (resp. paused or completed) or in "Reset Required" state (usually after an error).

Client Implementation

Because the server provides several endpoints, there are many ways to implement a client. We do not provide any client applications which are intended for real usage, however there are several demo clients that show how a client application can use the API. You can download the code of these demo clients at <https://github.com/BrukerSPR/Bruker-SPR-Remote-API>.

C#

The C# implementation can be found in the RemoteAppTestClient project. It is a Windows Forms application where the client object is based on the [System.ServiceModel.ClientBase](#) class using a [IClientChannel](#) that implements the same data contract interface as the server. The application provides controls to test almost all methods of the API and connects to the service using webHttp or TCP. The connection can be configured in the application's config file.

C++

Two options are given to create a C++ client. Option number 1 depends on Microsoft and its WebServices.dll. This way is described in the ReadMe.txt of the *RemoteAppTestClientConsoleNetTCP* project. The 2nd option is platform independent, uses [gSOAP](#) and is described in the ReadMe.txt of the *RemoteAppTestClientConsole* project. Both options do not have a hard-coded connection and their classes (*.c/*.cpp and *.h files) are generated with the described tools (see the specific ReadMe.txt).

Both demo applications are simple console applications providing the same features. Their usage is self-explanatory. After launching them they ask the user to enter a connection URI (or choose the default) and then offer different options which the user can select by entering the corresponding key.

Python, Java Script or other Web Languages

The webHttp endpoint can also be accessed using an arbitrary web browser. The data exchanged is based on JSON and requests can be submitted by just typing into the browser's address bar. It is hence easy to write a client in "web languages" like Python or Java Script. A python example can be found in the *RemoteAppTextClientPython* project. It contains the class *Mass1Api* which allows to call the remote API functions and a *Tools* class providing some helper methods. These two classes make it trivial to write a simple script to enqueue operations. An example for such a script is "RemoteAppTestClientPython.py".

Interface Description

To keep the client implementation as flexible as possible, our server does not use callbacks. This means that the server is not capable of contacting the client framework actively if an error has occurred or the state of the software has changed. The client must ask at regular intervals instead (polling). Messages by the control software are queued in a thread-safe queue instead. The arguments and return types are restricted to very fundamental types like numbers, Boolean values, strings, and arrays.

Methods

In Table 1 all available API methods are listed and described.

Remark: Whenever “method” (or “methods”) is part of the name of the API method or its description, not the API method itself is meant but the methods that can be created and run with the Control Software (i.e. are part of the runset of the Control Software). These methods may be also named procedures.

Table 1: API methods

Method Name	Return Type	Arguments	Description
Get information about methods			
GetNamesOfMethods	string[]	-	Returns an array of all method names that are loaded in the current runset library.
GetNamesOfMethodsOfAssayType	string[]	string assayType	Returns an array containing the names of all methods of the provided assay type that are loaded in the current runset library.
GetAssayTypesOfAllMethods	string[]	-	Returns an array containing the assay types of all methods that are loaded in the current runset library.
GetNameOfCurrentMethod	string	-	Returns the name of the method that is currently running or selected for running. Returns an empty string if there is no method selected.
GetAssayTypeOfCurrentMethod	string	-	Returns the assay type of the method that is currently running or selected for running. Returns an empty string if there is no method selected.

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

GetAssayTypeOfMethod	string	string methodName	Returns the assay type of the method with the given name. Returns an empty string if there is no method with the given name in the runset library.
Get information about runsets			
GetNamesOfRunsets	string[]	-	Returns an array of all runset names that are loaded in the current runset library.
GetNamesOfRunsetsOfAssayType	string[]	string assayType	Returns an array containing the names of all methods of the provided assay type that are loaded in the current runset library.
GetAssayTypesOfAllRunsets	string[]	-	Returns an array containing the assay types of all runsets that are loaded in the current runset library.
GetNameOfCurrentRunset	string	-	Returns the name of the runset that is currently running or selected for running. Returns an empty string if there is no runset selected.
GetAssayTypeOfCurrentRunset	string	-	Returns the assay type of the runset that is currently running or selected for running. Returns an empty string if there is no runset selected.
GetAssayTypeOfRunset	string	string runsetName	Returns the assay type of the runset with the given name. Returns an empty string if there is no runset with the given name in the runset library.
GetMethodNamesOfRunset	string[]	string runsetName	Returns the names of the methods in the runset with the given name.
Prepare and run assays			

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

SelectMethod	bool	string methodName	Selects the specified method for running. The method is wrapped into a runset that can then be started using the StartSelectedRunset method. Returns 'false' if the method could not be selected. Check for error messages in that case. Possible reasons are that the operation mode doesn't allow to select a method or that the method with the provided name doesn't exist in the runset library.
SelectRunset	bool	string runsetName	Select the specified runset for running. The selected runset can then be started using the StartSelectedRunset method. Returns 'false' if the runset could not be selected. Check for error messages in that case. Possible reasons are that the operation mode doesn't allow to select a runset or that the runset with the provided name doesn't exist in the runset library.
CreateRunset	bool	string[] methodNames	Creates a runset that contains the methods specified by the passed array of method names and selects the runset for running so it can be started using the StartSelectedRunset method. Returns 'false' if the runset couldn't be created. Check for error messages in that case. Possible reasons are that the operation mode doesn't allow to select a runset or that one or more methods with the provided names do not exist in the runset library.
SetSamplePlateId	bool	int methodName, string plateId	This method can only be used with SPR-24pro and SPR-32pro and will throw a NotSupportedException otherwise! Assigns an ID to the exchangeable plate of the method with the specified index in the currently selected runset. Returns 'false' if method index is out of range or the method with the given index doesn't use a sample plate. Check for errors in this case.

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

GetSamplePlateId	string	int methodIndex	This method can only be used with SPR-24pro and SPR-32pro and will throw a NotSupportedException otherwise! Gets the ID of the sample plate of the method with the given index in the current runset. Returns an empty string if the index is out of range, if there is no runset, no plate, or the plate doesn't have an ID assigned.
GetCurrentSamplePlateId	string	-	This method can only be used with SPR-24pro and SPR-32pro and will throw a NotSupportedException otherwise! Gets the ID of the sample plate that is currently inside the machine. Returns an empty string if there is no plate or ID is unknown.
SetPlateId	bool	int methodIndex, string plateId, string deckLocation	Assigns an ID to the plate at the provided deck location of the method with the specified index in the currently selected runset. Throws an exception if the deck location does not exist. In order to avoid this, make sure that the provided deck location comes from the array that you receive using the GetDeckLocations method. Returns 'false' if method index is out of range or the method with the given index doesn't use a sample plate. Check for errors in this case.
GetPlateId	string	int methodIndex, string deckLocation	Gets the ID of the plate at the provided deck location of the method with the given index in the current runset. Throws an exception if the deck location does not exist. To avoid this, make sure that the provided deck location comes from the array that you receive using the GetDeckLocations method. Returns an empty string if the index is out of range, if there is no runset, no plate, or the plate doesn't have an ID assigned.

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

GetCurrentPlateId	string	string deckLocation	Get the ID of the plate at the provided deck location that is currently inside the machine. Throws an exception if the deck location does not exist. To avoid this, make sure that the provided deck location comes from the array that you receive using the GetDeckLocations method. Returns an empty string if there is no plate or ID is unknown.
GetDeckLocations	string[]	-	Gets the names of all deck locations of the machine: R, S, V for SPR-24pro and SPR-32pro, W, T, P1, P2, P3, P4 for SPR #64.
MoveSamplePlateTrayOut	bool	-	This method can only be used with SPR-24pro and SPR-32pro devices and will throw a <code>NotSupportedException</code> otherwise. Moves the sample plate tray out of the machine so that the plate mover can get or put a plate from it. Returns 'false' if this is not possible. This might be the case if there is a runset or maintenance procedure currently running or if a reset is required. Check for errors in this case.
MoveSamplePlateTrayIn	bool	-	This method can only be used with SPR-24pro and SPR-32pro devices and will throw a <code>NotSupportedException</code> otherwise. Moves the sample plate tray into the housing. Returns 'false' if this is not possible. This might be the case if there is a runset or a maintenance procedure is currently running or if a reset is required. Check for errors in this case.

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

MoveReagentPlateTrayOut	bool	-	This method can only be used with SPR-24pro and SPR-32pro devices and will throw a <code>NotSupportedException</code> otherwise. Moves the plate tray out of the machine so that the plate mover can get or put a plate from it. Returns 'false' if this is not possible. This might be the case if there is a runset or maintenance procedure currently running or if a reset is required. Check for errors in this case.
MoveReagentPlateTrayIn	bool	-	This method can only be used with SPR-24pro and SPR-32pro devices and will throw a <code>NotSupportedException</code> otherwise. Moves the plate tray into the housing. Returns 'false' if this is not possible. This might be the case if there is a runset or a maintenance procedure is currently running or if a reset is required. Check for errors in this case.
OpenSideDoor	bool	-	This method can only be used with SPR #64 device and will throw a <code>NotSupportedException</code> otherwise. Opens the side door at the right side of the device. Returns 'false' if this is not possible. This might be the case if there is a runset or maintenance procedure is currently running or if a reset is required. Check for errors in this case.
CloseSideDoor	bool	-	This method can only be used with SPR #64 device and will throw a <code>NotSupportedException</code> otherwise. Close the side door at the right side of the device. Returns 'false' if this is not possible. This might be the case if there is a runset or maintenance procedure is currently running or if a reset is required. Check for errors in this case.

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

StartSelectedRunset	bool	-	Starts the currently selected runset. Returns 'false' if this is not possible. Check for errors in this case. Possible reasons might be that the system is not ready, has an error, or no selected runset is selected.
StartSelectedRunsetFrom	bool	int methodIndex	Starts the selected runset starting from the method with the given index. Returns 'false' if this is not possible. Check for errors in this case. Possible reasons might be that the system is not ready, has an error, the index of the method is out of range, or no runset is selected.
PauseRunsetAfter	bool	int pauseMode	Pause the runset that is currently running. The pauseMode argument sets when the runset is paused (1 = after current command, 2 = after current cycle, 3 = after current method). Returns 'false' if this is not possible because there is currently no runset running.
ResumeRunset	bool	-	Resumes the paused runset where it was paused. Returns 'false' if this is not possible. Check for errors in this case. Reasons might be that the operation mode is different from 'paused' or the device is not ready or has an error.
ResetRunset	bool	-	Resets the status of the runset after it was interrupted. Returns 'false' if this is not possible because the runset is currently running or if there is no runset to reset. Running the reset method multiple times on the same runset is possible.

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

AbortScript	bool	-	Aborts the currently running script immediately. That can either be a runset or a maintenance command. Returns 'false' if this is not possible because there is nothing running or the running script is already a cancel script (since the device runs a script to clean up the flow cell when interrupted in the middle of an injection).
LeaveStandby	bool	-	Leaves the standby mode. Returns 'false' if the machine is not in standby mode when this API method is called.
SetStandbyAfterFinish	bool	bool goToStandby	Sets whether the device should go into standby mode after the runset has finished. Always return true.
GetStandbyAfterFinish	bool	-	Gets whether the device should go into standby mode after the runset has finished. Returns 'true' = yes or 'false' = no.
Maintenance			
GetNamesOfMaintenanceProcedures	string[]	-	Returns an array containing the names of all maintenance commands that are available and can be run without user interaction.
RunMaintenanceProcedure	bool	string procedureName	Starts the maintenance procedure with the given name. Returns 'false' if this is not possible. Check for errors in this case. Reasons might be that the device is not ready or there is already something running.
Status			
GetOperationMode	int	-	Gets the current operation mode of the device. The meaning of the return value can be found in Table 2.
IsChipDocked	int	-	Checks whether a chip is docked. Returns '1' = docked, '0' = undocked, '-1' = unknown.

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

IsSamplePlateTrayIn	int	-	Check whether the sample plate tray is inside the machine. Returns '1' = inside, '0' = outside, '-1' = unknown.
IsReagentPlateTrayIn	int	-	Check whether the reagent plate tray is inside the machine. Returns '1' = inside, '0' = outside, '-1' = unknown.
IsSideDoorOpen	int	-	This method can only be used with SPR #64 device and will throw a <code>NotSupportedException</code> otherwise. Checks whether the door at the right side of the device is open. Meaning of the return value is: 1 = open, 0 = closed, -1 = unknown.
HasMessage	bool	-	Checks whether the server's message loop contains at least one new message.
GetMessage	string[]	-	Gets the first message from the server's message queue. A message consists of three strings in an array where the first string is a time stamp formatted like yyyyMMdd-hhmmss.fff, the second gives information about the type of message (Error, Warning, DataSaved, OperationModeChanged, ...), and the third string is the message itself. Call this method on a regular base using polling and check the result for null or use the <code>HasNewMessage</code> method first. A complete list of error types can be found in Table 3.
HasErrors	bool	-	Check whether there are any errors that are not covered by the operation mode.
GetErrors	string[]	-	Gets an array of strings with information about errors like 'no data', 'runset error', ...
HasWarnings	bool	-	Check whether there are any warnings that are not covered by the operation mode.

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

GetWarnings	string[]	-	Gets an array of strings with information about warnings like 'temperature unstable', 'runset warning', ...
-------------	----------	---	---

Operation Modes

Table 2: Operation modes

Numeric Value	Meaning	Description
0	Idle	The device is ready to run a script which can be a runset or a maintenance procedure.
1	Paused	The device is idle but there is a runset selected that was paused during run and might be resumed. The ResetRunset API method must be called before a new runset can be started. Running a maintenance command in between is possible, however.
2	Completed	The device is idle and the selected runset has been completed. The ResetRunset API method must be called before a new runset can be started.
3	Interrupted	The device is idle but the selected runset was interrupted by an emergency script. The ResetRunset API method must be called before a new runset can be started.
10	StandBy	The device is currently in standby mode. You must leave standby to get into the idle, paused, or completed mode.
11	ResetRequired	This state occurs after a stepper motor had an error or if a script was aborted due to an error or by user interaction. The maintenance method 'Reset System' must be run with the API method RunMaintenanceProcedure to get into idle, paused, or completed mode.
12	MaintenanceRequired	The device is not ready for operation but needs technical maintenance. Currently not used.
13	DoorOpen	The front door of the device is open. It must be closed before the machine can be used.
20	Running	There is currently a script running. This can be either a runset or a maintenance method.

Error Types

Table 3: Error types

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

Type	Description
Error	A critical error has occurred; current operation cannot proceed
Warning	An uncritical error occurred
StateChanged	The operation mode has changed (see Table 2)
DataCleared	The recorded data were cleared from the application's memory
DataSaved	The recorded data were saved to file
RunsetStarted	A runset was started
RunsetCompleted	A runset has finished
MethodStarted	A method was started
MethodCompleted	A method has finished
CycleStarted	A cycle was started
CycleCompleted	A cycle has finished
InjectionStarted	An injection was started
InjectionCompleted	An injection has finished
Rollback	Events have been removed. This happens usually if a batch job was interrupted during an injection. In this case, the injection events are removed, and the injection is restarted when the job is resumed.
AppMessage	Other message from script

Device Geometry SPR-24pro and SPR-32pro

The geometry of the device is shown below. All distances are given in millimeters. The plate tray for the sample plate can be moved out to the right side of the device, so the plate can be exchanged by a plate robot controlled by a lab automation framework. The plate mover must be placed at the right side of the device. The space at the front of the machine should stay accessible for manual operations and maintenance work. There is also space for a plate on the left side of the device, however, the current software does not support to exchange it at runtime.

It is important to mark the position of the machine and the plate mover so their relative position can be restored after one of them has been moved during service work.

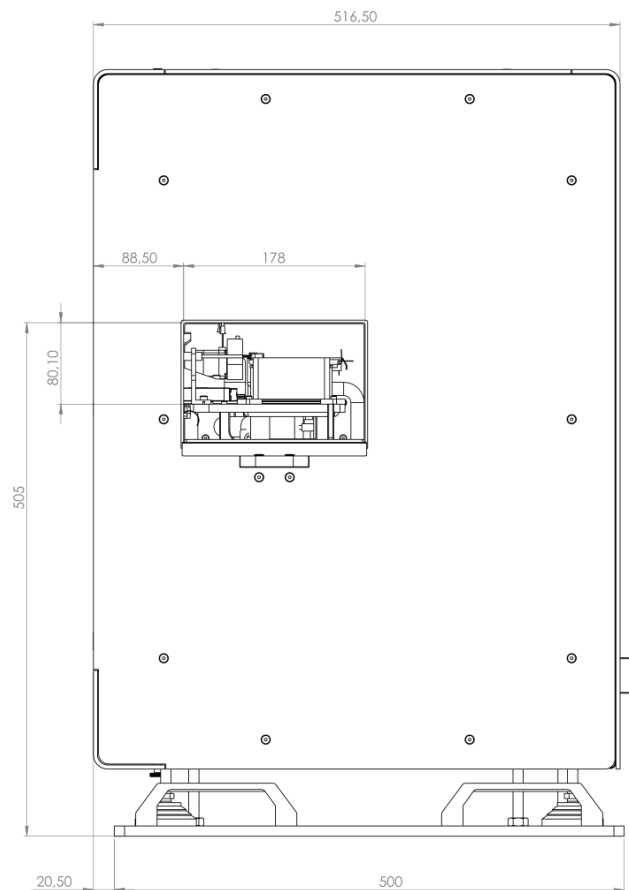


Figure 1: side view

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

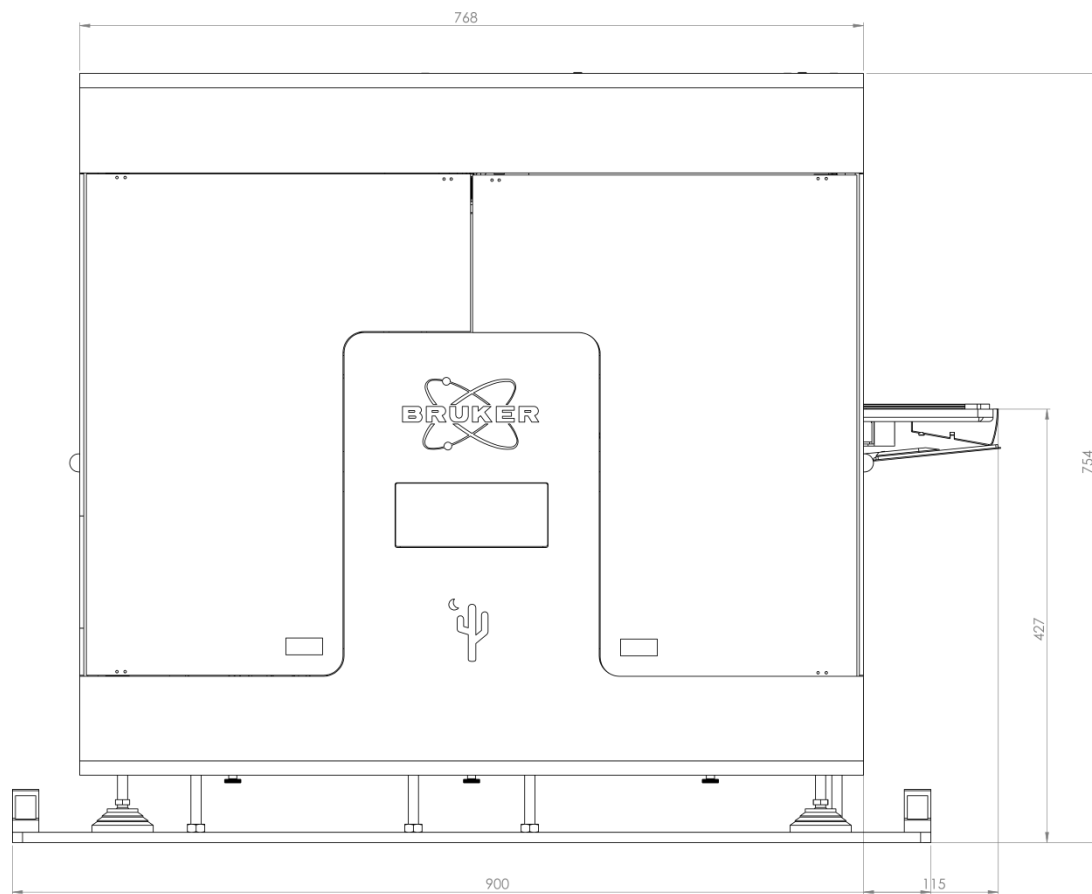


Figure 2: front view

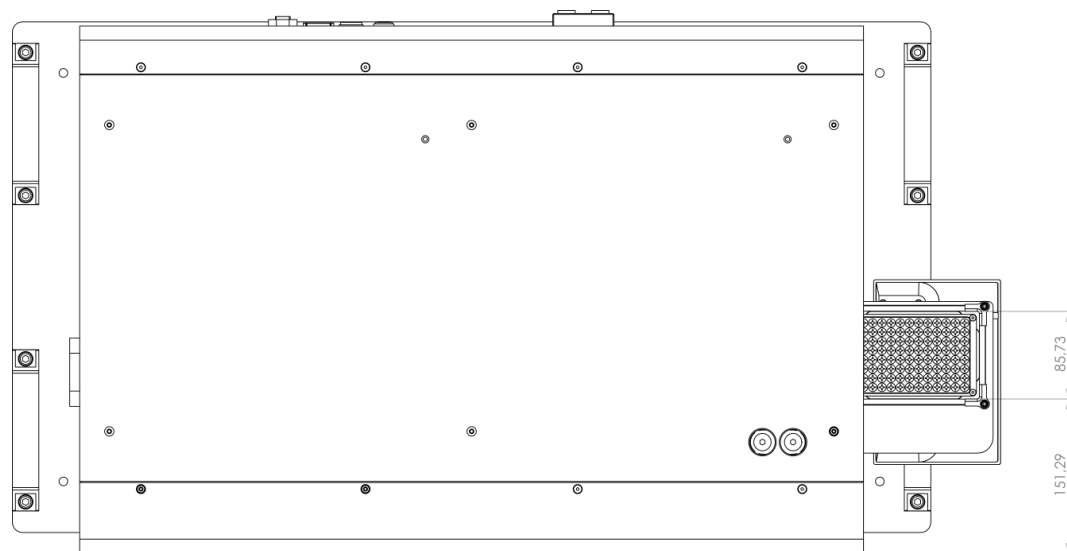


Figure 3: top view

Device Geometry SPR #64

The geometry of the device is shown below. All distances are given in millimeters. The plate mover must be placed at the right side of the device. The space at the front of the machine should stay accessible for manual operations and maintenance work.

It is important to mark the position of the machine and the plate mover so their relative position can be restored after one of them has been moved during service work.

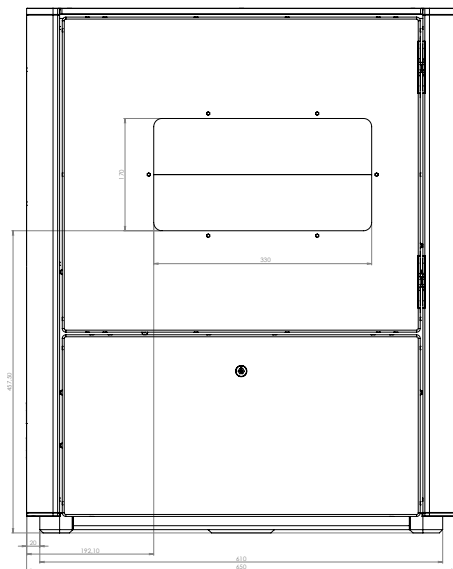


Figure 4: side view

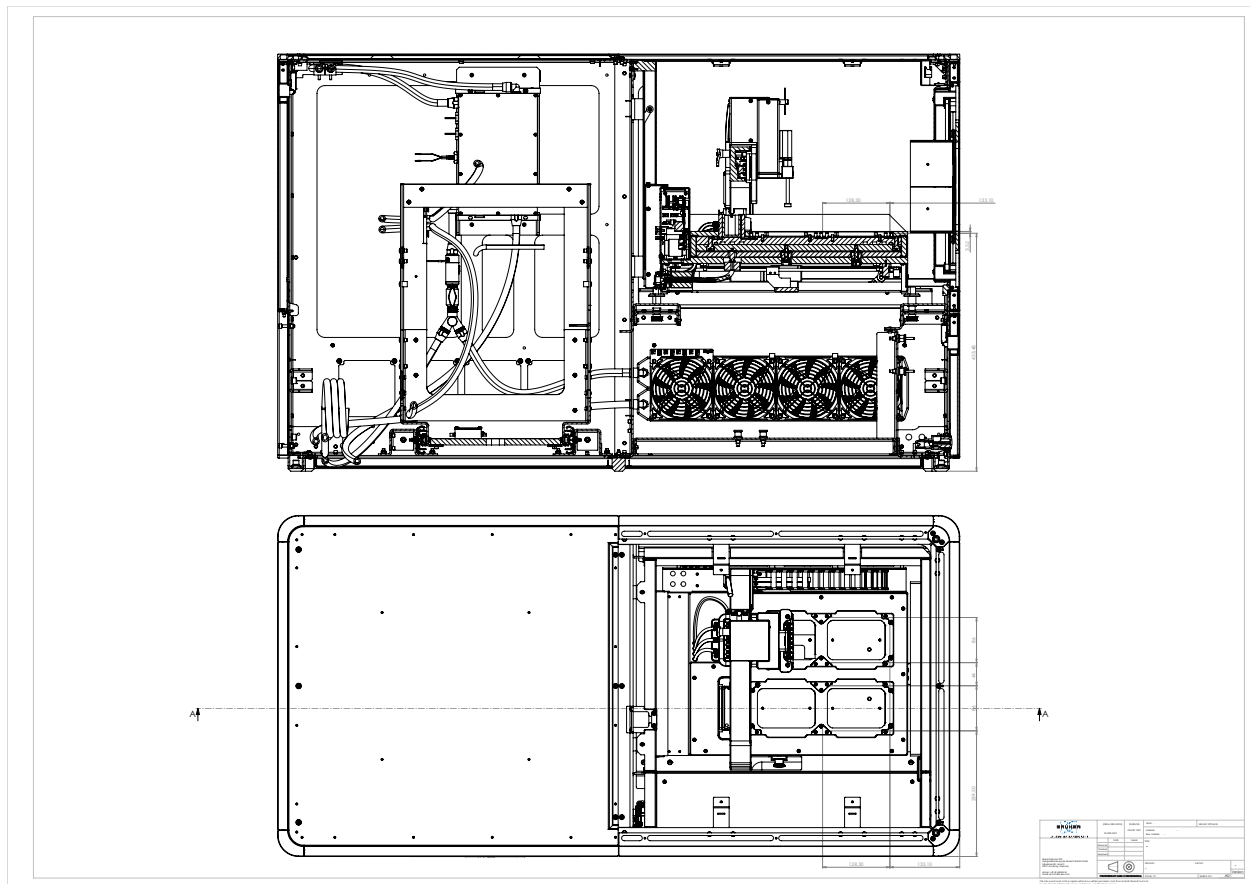


Figure 5: front view and top view

Bruker Biosensors

Bruker Daltonics SPR, Falkenried 88, 20251 Hamburg, Germany
Zweigniederlassung der Bruker Daltonics GmbH & Co. KG

Manufacturer

Manufactured by:



Bruker Daltonics GmbH & Co KG
Fahrenheitstr. 4
28359 Bremen
Germany
Phone: +49 (421) 2205-0

Revision Notes

Revision ID	Date	Changes	Name
A	2019-11-04	Initial revision for version 3.5.1	Dirk Dammann
B	2019-11-29	Replaced images	Dirk Dammann
C	2020-06-05	Added ini configuration	Dirk Dammann
D	2024-04-22	Extended for version 3.9 supporting SPR #64	Dirk Dammann
E	2024-07-24	Added methods about moving the reagent plate in and out	Dirk Dammann
F	2026-01-09	Replaced SPR #64 images and some minor clarifications	Dirk Dammann
G	2026-01-09	Updated header	Sven Malik
H	2026-08-18	Clarified ini settings and fixed wrong links in Endpoint Configuration.	Tobias Ernst